

SonoTAP® II


*Echogenic, Non-Stimulating
Single Shot Plane Block Needle*

Plane blocks like TAP, QL, ESP etc. are gaining momentum and are a valuable alternative to Epidural Anesthesia or General Anesthesia. However, these blocks rely solely on users experience and needle shaft and tip echogenicity. With SonoTAP II, PAJUNK® has developed a special plane block needle for ultrasound guidance only, which includes the patented Cornerstone Reflectors for optimized needle echogenicity.




Non-Removable Injection Tubing

- Closed system
- No leakage

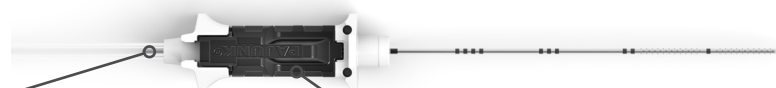


Colour-Coded Needle Hub

- Easy identification of needle size

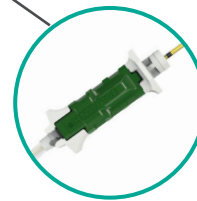



Depth Markings

- Easy to read and identify

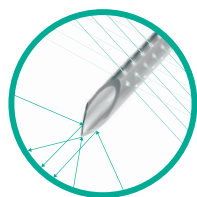

Echogenic Needle Tip

- Facet tip with two inclination angles
- Improves needle tip visibility under ultrasound

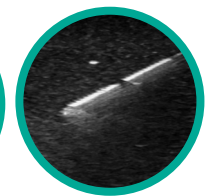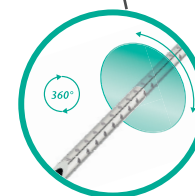

Cornerstone Reflectors

- 360 degree graduations on the first 20mm of the needle.
- Optimised ultrasound visibility of needle shaft¹
- Reliable and optimised needle echogenicity at higher angles

1. Edgcombe H., Hocking G. (2010 March-April) Sonographic identification of needle tip by specialists and novices: a blinded comparison of 5 regional block needles in fresh human cadavers. Reg. Anesth. Pain Med., 35(2): 207-211

SonoTAP® II

Echogenic, Non-Stimulating
Single Shot Plane Block Needle

SonoTAP II	Größe	Artikel-Nr.	NRFit® Artikel-Nr.	VE
Facettenschliff	24G x 40mm (1 5/8")	1285-3Y040	1265-3Y040	10
	22G x 50mm (2")	1285-3E050	1265-3E050	10
	22G x 80mm (3 1/5")	1285-3E080	1265-3E080	10
	21G x 110mm (4")	1285-3F110	1265-3F110	10
	21G x 150mm (6")		1265-3F150	10

NRFit® is a registered trademark of GEDSA used with their permission.

PAJUNK® GmbH Medizintechnologie
Global Headquarter
Karl-Hall-Straße 1
78187 Geisingen, Germany
Tel.: +49 7704 9291-0
Fax: +49 7704 9291-600
info@pajunk.com
pajunk.com

PAJUNK® Medical Produkte GmbH
Central Europe
Karl-Hall-Straße 1
78187 Geisingen, Germany
Tel.: +49 7704 8008-0
Fax: +49 7704 8008-150
info@pajunk-medipro.com
pajunk.de & pajunk.eu

PAJUNK® Swiss AG
Switzerland
Neuhofstrasse 5A
6340 Baar, Switzerland
Tel.: +41 41 760 96 27
info@pajunk.ch
pajunk.ch

PAJUNK® Medical System L.P.
United States of America
4575 Marconi Drive
Alpharetta, GA 30005, USA
Tel.: +1 770.493.6832
Fax: +1 678.514.3388
info@pajunk-usa.com
pajunkusa.com

PAJUNK® UK Medical Products Ltd
United Kingdom
Unit D1 The Waterfront
Goldcrest Way, Newburn Riverside
Newcastle upon Tyne, NE15 8NY, UK
Tel.: +44 191 264 7333
info@pajunk.co.uk
pajunk.co.uk